



Baropodometric comparison of orthopedic footwear to assess the effectiveness of pairs of orthopedic shoes for reducing the forefoot pressure

Hasan Basri Sezer, MD¹, Alexandre Hardy, MD¹, Ronny Lopes, MD²

¹Department of Sports Surgery, Clinique Du Sport Paris V, Chirurgie Du Sport, Paris, France

²Department of Sports Surgery, Health Atlantic, Chirurgie Du Sport, Saint Herblain, France

Orthopedic shoes are mainly used for the conservative treatment of diabetic foot and following injuries or foot surgery. The goal of using orthopedic footwear is to reduce plantar pressure, particularly on the forefoot. Historically, this has involved a single shoe with the forefoot weight offloaded by the heel, with or without the support of the forefoot.^[1] A more modern design has recently been marketed with pairs of offloading shoes with full, convex, rigid soles, in other words, rocker bottom soled shoes.

Foot surgery is in continuous evolution. Progress in surgical techniques of the forefoot has significantly changed the practices and goals of orthopedic footwear.^[2,3] Greater understanding and the more extensive use of percutaneous or minimally invasive surgery associated with increasingly reliable

Received: May 04, 2022

Accepted: June 07, 2022

Published online: July 06, 2022

Correspondence: Hasan Basri Sezer, MD. Clinique Du Sport Paris, Chirurgie Du Sport, 75005 Paris, France.

E-mail: mrhbs@hotmail.com

Doi: 10.52312/jdrs.2022.702

Citation: Sezer HB, Hardy A, Lopes R. Baropodometric comparison of orthopedic footwear to assess the effectiveness of pairs of orthopedic shoes for reducing the forefoot pressure. Jt Dis Relat Surg 2022;33(2):323-329.

©2022 All right reserved by the Turkish Joint Diseases Foundation

This is an open access article under the terms of the Creative Commons Attribution-NonCommercial License, which permits use, distribution and reproduction in any medium, provided the original work is properly cited and is not used for commercial purposes (<http://creativecommons.org/licenses/by-nc/4.0/>).

ABSTRACT

Objectives: This study aimed to evaluate the distribution of foot pressure with the most frequently used orthopedic shoes and demonstrate the effect of offloading philosophy on the pressure distributions of rocker bottom or heel support shoes applied unilaterally or bilaterally.

Materials and methods: Three bilateral and four unilateral, a total of seven shoe designs with sensors included in the insole were tested by the same subject in a standard acquisition protocol. Two of the unilateral and one of the bilateral shoes had heel support, while others had a rocker bottom design. A descriptive analysis was performed for each shoe and compared to a reference value obtained from a standard shoe.

Results: Shoes with an offloading heel resulted in a greater reduction in pressure on the forefoot than other models. All other shoes increased pressure on the first metatarsophalangeal joint. Heel offloading models performed the best in forefoot offloading, and bilateral heel offloading shoes performed the best in first metatarsal offloading, with the highest scores of 83% and 82%, respectively.

Conclusion: This study showed that orthopedic shoes sold in pairs could reduce pressure on the forefoot at a comparable level to unilateral shoes. It supports their use to limit the disadvantages of single orthopedic shoes, such as limb length discrepancies.

Keywords: Forefoot relief, offloading, plantar pressure distribution, postoperative shoe.

osteosynthesis techniques have resulted in earlier postoperative weight bearing and early discharge from the hospital.^[4,5] Although there exist numerous studies that evaluate the distribution of plantar pressure with different orthopedic shoes^[1,6-13] and its role in the diabetic foot,^[14,15] to our knowledge, none of them compare pairs of orthopedic footwear to other single orthopedic shoes. Sold in pairs, these shoes have the capacity to limit the limb length differences. Overcoming the leg length discrepancy

might increase the comfort and speed of walking while decreasing the postoperative falling rates.^[16-18] However, prescription of bilateral shoes is not a part of daily practice, probably due to economic concerns and refusal of refund by the assurance for the healthy side.

In 2019, Dearden et al.^[19] performed a prospective study and did not find any difference in the clinical and radiological progress following surgery of the forefoot between heel offloading shoes and rocker bottom shoes. However, they reported greater comfort and patient satisfaction and better stability with full rocker bottom soled shoes. A reasonable explanation may be the theoretical advantage of bilateral shoes, which may help patients adapt to a relatively new postsurgical condition.

The main goal of this study was to evaluate the distribution of foot pressure with the most frequently used orthopedic shoes and demonstrate the effect of offloading philosophy on the pressure distributions of rocker bottom or heel support shoes applied unilaterally or bilaterally. With the use of unilateral postoperative shoes, there is presumed to be a difference between the two sides in terms of elasticity of footwear, comfort, and walking patterns, although the patients wear a suitable shoe on the healthy side. We hypothesized that the reduction in plantar pressure with orthopedic footwear sold in pairs would be superior to one-sided shoes that are commercially available.

MATERIALS AND METHODS

A protocol was designed for data collection. The left foot of a 25-year-old female with no past medical

or surgical history was bandaged as if she had undergone forefoot surgery. The subject was 1.64 meters tall, weighed 52 kg, and wore a size 38 shoe. The reference shoe was a comfortable, standard sports shoe. When unilateral orthopedic footwear was being tested on the left foot, the reference shoe was worn on the right foot without a sock.

The subject walked 10 m in a straight line in each tested shoe; this was performed twice. The first 10 m allowed the subject to become familiar with walking in the shoe, and the second was used to register the data. Thus, the second walk constituted the data. Walking speed was natural and unguided. Recording of plantar pressure patterns was obtained with an insole (Moticon Science, Munich, Germany; Figure 1). This 3 mm insole was bilaterally placed. Walking was recorded, and an automatic video synchronization program associated the measurements of plantar pressure with the gait cycle. The amount of pressure exerted was quantitatively measured (N/cm^2) and retranscribed to color code. The same insoles were used in all of the tested shoes.

Seven shoe brands were compared in this study. Podalux™ (DJO Global, Lewisville, TX, USA), CHV (Sober, Saint-Georges-de-Reneins, France), Barouk® second-generation (Romans Industrie, Romans, France), and Orthowedge (DARCO (Europe) GmbH, Raisting, Germany) shoes were tested unpaired. Halten (Podonov, Saulx-les-Chartreux, France), Gemini, and Gemini Alpha® (Implant Service Orthopedie, Ris-Orangis, France) shoes were tested in pairs. Barouk, Orthowedge, and Gemini Alpha shoes have heel support, while others have rocker bottom soles. The characteristics of the seven tested shoes are presented in Table I and Figure 2.

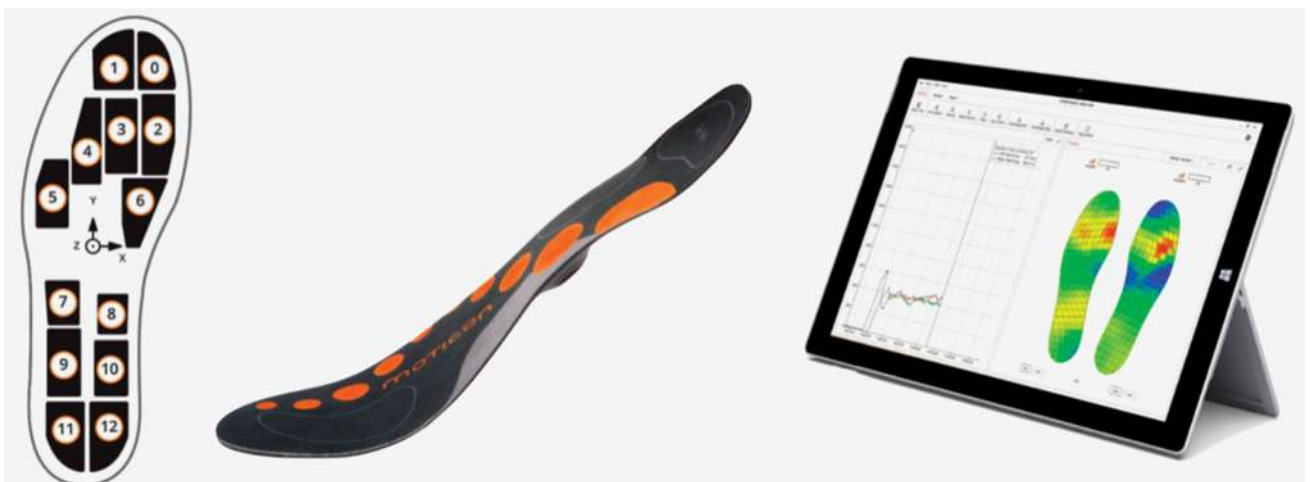


FIGURE 1. Recording system of plantar pressure.

TABLE I Characteristics of tested shoes					
Name	Number of shoes	Type of weight-bearing	Waterproof	Anterior protection	Fastener
Podalux (DJO Global)	1	Rocker bottom sole	No	Yes	Velcro
CHV (SOBER)		Rocker bottom sole	No	Yes	Velcro
Barouk (ROMANS)		Heel	Yes	No	Velcro
Orthowedge (DARCO)		Heel	No	No	Velcro
Halten (PODONOV)	2	Rocker bottom sole	Yes	No	Laces
Gemini (ISO)		Rocker bottom sole	Yes	Yes	Laces + double zippers
Gemini Alpha (ISO)		Heel	Yes	Yes	Laces + double zippers

Statistical analysis

All analyses were performed using the IBM SPSS version 23.0 software (IBM Corp., Armonk, NY, USA). A descriptive analysis of the measured plantar pressure values was performed. The number of steps, walking speed in terms of strides per minute, mean forefoot pressures (N/cm²), and mean first metatarsophalangeal (MTP1) pressures (N/cm²) were

measured and compared. The insole consisted of seven forefeet and six hindfeet, a total of 13 pressure sensors, and a three-dimensional acceleration sensor with an embedded ANT radio source, a thermal sensor, a flash memory, and a power supply. The S2 and S3 sensors represented the contact area of the MTP1, and the mean values of these sensors were utilized to express the corresponding parameter.

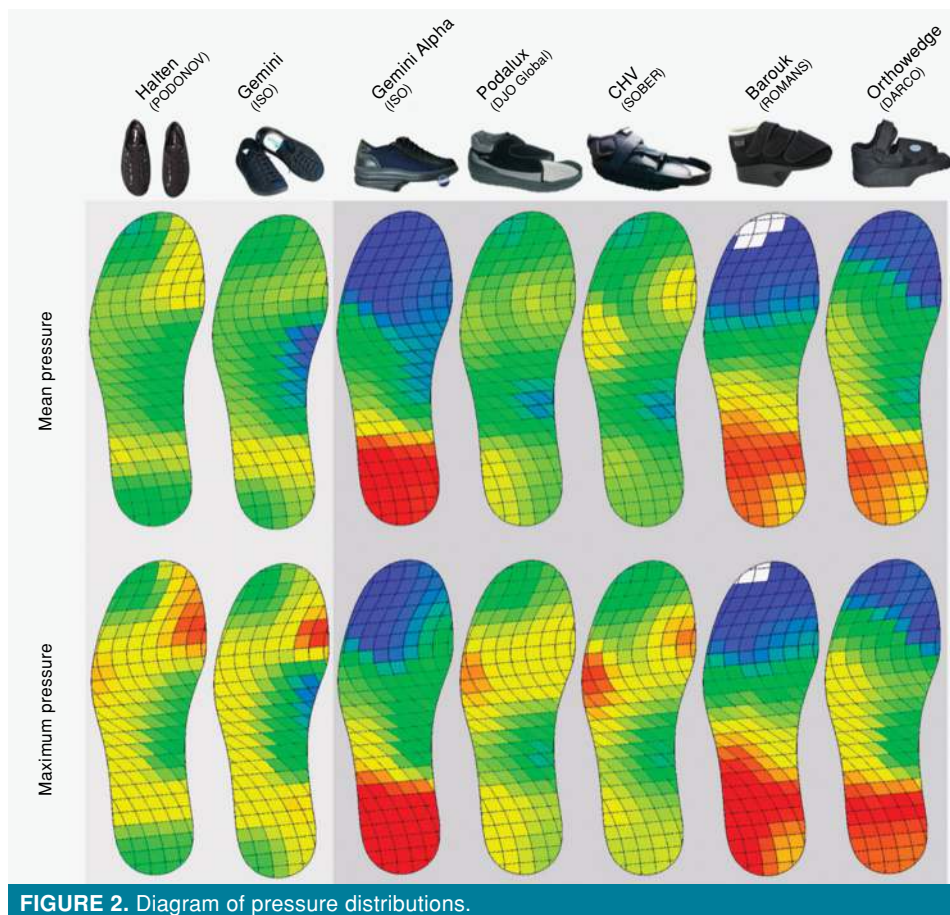


FIGURE 2. Diagram of pressure distributions.

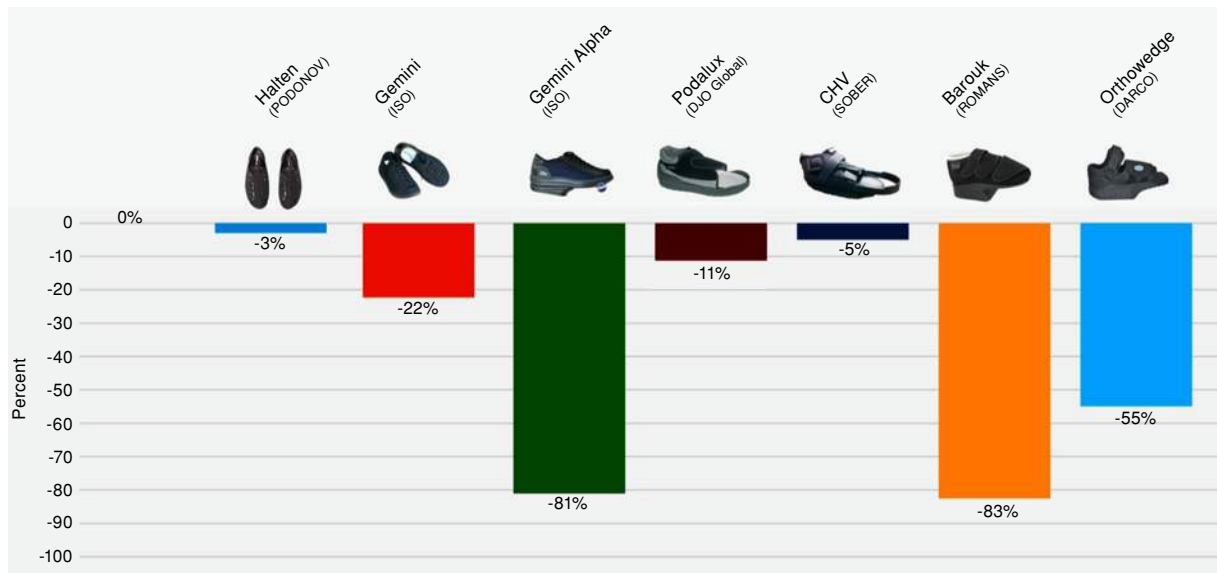


FIGURE 3. Percentage of offloading of the forefoot according to shoe design.

The insole was capable of wireless transfer with the capacity of providing a visual data interface (Figures 1 and 2).

RESULTS

All measured values of the plantar pressures are presented in Table II and Figure 2. The paired orthopedic footwear performed a similar reduction in plantar pressure to unilateral orthopedic shoes. This was against our initial hypothesis. Walking speed was similar in the tested shoes. Walking was slower in shoes with a forefoot offloading heel or a single shoe. Shoes with an offloading heel, namely Barouk second-generation, Gemini Alpha, and the Orthowedge,

resulted in the greatest reduction in pressure of the forefoot compared to the reference shoe with 83%, 81%, and 55%, respectively. Among the shoes with a full rocker bottom sole, the Gemini shoe reduced pressure on the forefoot by 22%, which surpassed the values produced by Podalux, CHV, and Halten shoes. Figure 3 presents the reduction in forefoot pressure with each shoe.

If the mean pressure values of the S2 and S3 sensors were separately evaluated, representing the pressures on the MTP1, heel offloading property of Gemini Alpha, Barouk second-generation, and the Orthowedge shoes resulted in a reduction of 82%, 80%, and 73%, respectively, compared to the reference shoe.

Shoe	Steps	Speed (Strides/ minute)	Mean forefoot pressure (S0-S6)		Mean MTP1 pressure (S2+S3)	
			(N/cm ²)	%	(N/cm ²)	%
Reference shoe (ADIDAS)	45	53.55	21.26	100	5.83	100
Halten (PODONOV)	32	59.83	20.64	97	8.34	143
Gemini (ISO)	30	62.03	16.53	78	6.78	116
Gemini alpha (ISO)	40	54.58	4.03	19	1.07	18
Podalux (DJO Global)	32	59.28	18.88	89	6.65	114
CHV (SOBER)	33	59.00	20.20	95	7.65	131
Barouk (ROMANS)	45	53.55	3.72	17	1.19	20
Orthowedge (DARCO)	40	60.85	9.60	45	1.56	27

MTP1: Mean first metatarsophalangeal.

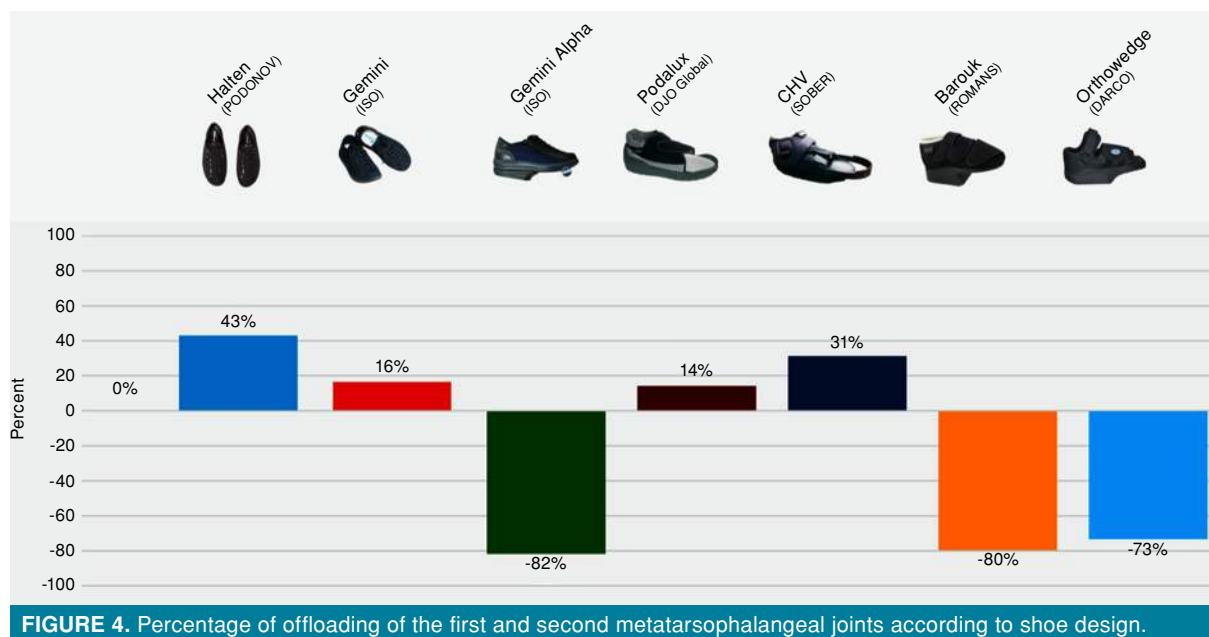


FIGURE 4. Percentage of offloading of the first and second metatarsophalangeal joints according to shoe design.

In addition, all of the shoes with full rocker bottom soles increased pressure on the MTP1 in different proportions. Gemini and Podalux shoes resulted in an increase of pressures by 16% and 14%, respectively. More remarkable increases in the MTP1 pressure were noted with CHV and Halten shoes by 31% and 43%, respectively. Figure 4 presents measurements of the pressure on the MTP1 with different shoes.

DISCUSSION

The knowledge on the characteristics of postoperative shoes continues to accumulate. As expected, orthopedic shoes sold in pairs may limit the problems of limb length discrepancies created by wearing a single orthopedic shoe.^[16-18,20,21] There is also existing evidence for the use of heeled shoes, which provide significantly more forefoot offloading (between 2.5 and 4 times more depending on the design) compared to rocker bottom soles.^[22] The findings of this study support the use of heel offloading shoes to acquire the highest degree of forefoot offloading, although rocker bottom shoes provided more speed in pairs.

Back and hip pain are the most frequently reported effects of the limb length discrepancies created by unequal orthopedic shoes.^[18] Michalik et al.^[16] compared a rocker bottom sole shoe to an offloading heel shoe and found a statistically significant difference in the oblique angle of the pelvis and the lateral deviation of the spine

compared to a regular shoe. Furthermore, they found hardly any difference in gait between the two orthopedic sole designs that were evaluated. They concluded that whatever the design of the shoe, a compensating shoe on the contralateral side should be used to transfer the pressure to the lateral side of the operated heel.^[16] Pairs of orthopedic shoes may not only be beneficial for proper transfer of load but also decrease falling risk in the elderly subjects caused by unbalanced shoes.^[23]

In this study, heeled shoes resulted in a remarkable reduction in the MTP1 pressure in a range of 73% to 82%. On the contrary, shoes with rocker bottom soles increased the forefoot pressure. Comparable results were reported by Van Schie et al.^[22] in 14 of the 17 patients evaluated by using the Darco® VFE (DARCO (Europe) GmbH, Raisting, Germany) shoe. Schuh et al.^[24] evaluated the characteristics of five shoes on the market and five prototypes to a reference shoe and found similar results. They blamed the increased torsional force on the hallux due to the excessive stiffness of the soles and assumed that a combination of excessive rigidity and an exceedingly posterior rocking line could have created this paradoxical effect. They concluded that the quality of orthopedic footwear was based on a subtle balance of the rigidity of the sole, the rolling characteristics, the rocking line, and shoe comfort. Moreover, they reported a need for a proper correlation between the length and the width (length-width ratio) to ensure patient comfort. There

were certain differences in the quantitative pressure values (N/cm²) reported in the study of Schuh et al.^[24] and the results of the present study. This difference was probably due to the system of measurement and calibration.^[25] The measurement system utilized in our study was validated by numerous studies.^[26-28] To our knowledge, this is the first study to evaluate the distribution of plantar pressures in orthopedic shoes sold in pairs.^[29]

The main limitation of this study is the evaluation of only a single healthy subject. The subject on whom the tests were performed was healthy and had no problem with her foot or limping. We did not report the measurements in the healthy foot, as this was not pertinent to our analysis. Second, the speed of walking was not guided, which has been reported to modify the plantar pressure.^[30] As the assessment of walking speed is not easy in a short distance, we assumed that the speed was nearly constant for 10 m and that any variations and their influence on the pressure measurements was negligible. Finally, postoperative bandaging following foot surgery markedly varies from one surgeon to the other and may significantly modify the distribution of the plantar pressure.

In conclusion, the reduction of pressure on the forefoot was significantly greater with a forefoot offloading heel support model compared to shoes with rocker bottom soles, which conversely increased the pressure on the first metatarsophalangeal joint. Paired heel supporting models may be considered as an alternative option during the early postoperative period to obtain significant forefoot offloading and comfort.

Ethics Committee Approval: The study protocol was approved by the IRB Elsan Ethical Committee (IRB#2021-07-LOPES-02).

Data Sharing Statement: The data that support the findings of this study are available from the corresponding author upon reasonable request.

Author Contributions: Writer, reviewer, analysis and interpretation: H.B.S.; Idea, concept, design of the study and supervisor: A.H.; Idea, concept, data collection and writer: R.L.

Conflict of Interest: The authors declared no conflicts of interest with respect to the authorship and/or publication of this article.

Funding: The authors received no financial support for the research and/or authorship of this article.

REFERENCES

- Deleu PA, Leemrijse T, Vandeleene B, Maldague P, Devos Bevernage B. Plantar pressure relief using a forefoot offloading shoe. *Foot Ankle Surg* 2010;16:178-82.
- Caravelli S, Mosca M, Massimi S, Costa GG, Lo Presti M, Fuiano M, et al. Percutaneous treatment of hallux valgus: What's the evidence? A systematic review. *Musculoskelet Surg* 2018;102:111-7.
- Maffulli N, Longo UG, Marinozzi A, Denaro V. Hallux valgus: Effectiveness and safety of minimally invasive surgery. A systematic review. *Br Med Bull* 2011;97:149-67.
- Robinson AH, Limbers JP. Modern concepts in the treatment of hallux valgus. *J Bone Joint Surg [Br]* 2005;87:1038-45.
- Mouton A, Le Strat V, Medevielle D, Kerroumi Y, Graff W; French association of foot surgery (AFCP). Patient's satisfaction after outpatient forefoot surgery: Study of 619 cases. *Orthop Traumatol Surg Res* 2015;101(6 Suppl):S217-20.
- Alexander IJ, Chao EY, Johnson KA. The assessment of dynamic foot-to-ground contact forces and plantar pressure distribution: A review of the evolution of current techniques and clinical applications. *Foot Ankle* 1990;11:152-67.
- Carl HD, Pfander D, Swoboda B. Assessment of plantar pressure in forefoot relief shoes of different designs. *Foot Ankle Int* 2006;27:117-20.
- Glod DJ, Fettinger P, Gibbons RW. A comparison of weightbearing pressures in various postoperative devices. *J Foot Ankle Surg* 1996;35:149-54.
- Lorei T, Klärner H, Rosenbaum D. Influence of postoperative shoes on plantar pressure patterns. *Z Orthop Ihre Grenzgeb* 2006;144:153-7.
- Nawoczinski DA, Birke JA, Coleman WC. Effect of rocker sole design on plantar forefoot pressures. *J Am Podiatr Med Assoc* 1988;78:455-60.
- Notni A, Fuhrmann RA. Measuring plantar weight distribution of shoes with forefoot support. *Z Orthop Ihre Grenzgeb* 1999;137:280-3.
- Schaff PS, Cavanagh PR. Shoes for the insensitive foot: The effect of a "rocker bottom" shoe modification on plantar pressure distribution. *Foot Ankle* 1990;11:129-40.
- Trnka HJ, Parks BG, Ivanic G, Chu IT, Easley ME, Schon LC, et al. Six first metatarsal shaft osteotomies: Mechanical and immobilization comparisons. *Clin Orthop Relat Res* 2000;3:256-65.
- Chapman JD, Preece S, Braunstein B, Höhne A, Nester CJ, Brueggemann P, et al. Effect of rocker shoe design features on forefoot plantar pressures in people with and without diabetes. *Clin Biomech (Bristol, Avon)* 2013;28:679-85.
- Bus SA, van Deursen RW, Armstrong DG, Lewis JE, Caravaggi CF, Cavanagh PR; International Working Group on the Diabetic Foot. Footwear and offloading interventions to prevent and heal foot ulcers and reduce plantar pressure in patients with diabetes: A systematic review. *Diabetes Metab Res Rev* 2016;32 Suppl 1:99-118.
- Michalik R, Siebers H, Claßen T, Gatz M, Rohof B, Eschweiler J, et al. Comparison of two different designs of forefoot off-loader shoes and their influence on gait and spinal posture. *Gait Posture* 2019;69:202-8.
- Meyer PE, Petersen D. Compensation in leg length inequality with orthopedic shoe measures. *Orthopade* 1992;21:174-83.
- Friberg O. Clinical symptoms and biomechanics of lumbar spine and hip joint in leg length inequality. *Spine (Phila Pa 1976)* 1983;8:643-51.
- Dearden PMC, Ray RI, Robinson PW, Varrall CR, Goff TJ, Fogarty KA, et al. Clinical and radiological outcomes of forefoot offloading versus rigid flat shoes in patients

- undergoing surgery of the first ray. *Foot Ankle Int* 2019;40:1189-94.
20. Paton JS, Thomason K, Trimble K, Metcalfe JE, Marsden J. Effect of a forefoot off-loading postoperative shoe on muscle activity, posture, and static balance. *J Am Podiatr Med Assoc* 2013;103:36-42.
 21. Wünnemann M, Klein D, Rosenbaum D. Effects of the Twin Shoe (Darco) to compensate height differences in normal gait. *Gait Posture* 2011;33:61-5.
 22. van Schie C, Ulbrecht JS, Becker MB, Cavanagh PR. Design criteria for rigid rocker shoes. *Foot Ankle Int* 2000;21:833-44.
 23. Hijmans JM, Geertzen JH, Dijkstra PU, Postema K. A systematic review of the effects of shoes and other ankle or foot appliances on balance in older people and people with peripheral nervous system disorders. *Gait Posture* 2007;25:316-23.
 24. Schuh R, Trnka HJ, Sabo A, Reichel M, Kristen KH. Biomechanics of postoperative shoes: Plantar pressure distribution, wearing characteristics and design criteria: A preliminary study. *Arch Orthop Trauma Surg* 2011;131:197-203.
 25. Murphy DF, Beynnon BD, Michelson JD, Vacek PM. Efficacy of plantar loading parameters during gait in terms of reliability, variability, effect of gender and relationship between contact area and plantar pressure. *Foot Ankle Int* 2005;26:171-9.
 26. Braun BJ, Veith NT, Hell R, Döbele S, Roland M, Rollmann M, et al. Validation and reliability testing of a new, fully integrated gait analysis insole. *J Foot Ankle Res* 2015;8:54.
 27. Braun BJ, Bushuven E, Hell R, Veith NT, Buschbaum J, Holstein JH, et al. A novel tool for continuous fracture aftercare - Clinical feasibility and first results of a new telemetric gait analysis insole. *Injury* 2016;47:490-4.
 28. Braun BJ, Pelz P, Veith NT, Rollmann M, Klein M, Herath SC, et al. Long-term pathological gait pattern changes after talus fractures - dynamic measurements with a new insole. *Int Orthop* 2018;42:1075-82.
 29. Atik OŞ. Which articles do the editors prefer to publish? *Jt Dis Relat Surg* 2022;33:1-2.
 30. Hennig EM, Staats A, Rosenbaum D. Plantar pressure distribution patterns of young school children in comparison to adults. *Foot Ankle Int* 1994;15:35-40.